

- Eimüller, G. J. Kusinski, "High resolution X-ray microscopy with magnetic contrast," *Journal de Physique IV*, in press.
7. P. Fischer, G. Schütz, G. Schmahl, P. Guttmann and D. Raasch, Imaging of Magnetic Domains with the X-ray microscope at BESSY using X-ray Magnetic Circular Dichroism, *Zf. Physik B* **101** (1996) 313–316.
8. P. Fischer, T. Eimüller, G. Schütz und G. Denbeaux, Imaging magnetic domain structures with Soft X-ray Microscopy, *J. Structural Chemistry* **14**(1) (2003) 39–47.
9. P. Fischer, T. Eimüller, G. Schütz, G. Denbeaux, A. Lucero, L. Johnson, D. Attwood, S. Tsunashima, M. Kumazawa, N. Takagi, M. Köhler, and G. Bayreuther, Element-specific imaging of magnetic domains at 25nm spatial resolution using soft X-ray microscopy, *Rev. Sci. Instr.* **72**(5) (2001) 2322–24.
10. G. Denbeaux, P. Fischer, G. Kusinski, M. Le Gros, A. Pearson, and D. Attwood, "A full field transmission X-ray microscope as a tool for high-resolution magnetic imaging," *IEEE Transactions on Magnetics*, **37**, 2764–2766, 2001.
11. P. Fischer, G. Denbeaux, T. Ono, T. Okuno, T. Eimüller, D. Goll, and G. Schütz, Study of magnetic domains with magnetic soft X-ray transmission microscopy, *J. Phys. D Appl. Phys.* **35**(19) (2002) 2391–97.
12. W. T. Lynch, L. A. Arledge, *Proc. Mater. Res. Soc. Symp.* **514** (1998) 11.
13. G. Schneider, G. Denbeaux, E. H. Anderson, B. Bates, A. Pearson, M. A. Meyer, E. Zschech and E. A. Stach, *Appl. Phys. Lett.* **81**, 14 (2002) 2535–2537.
14. G. Schneider, G. Denbeaux, E. Anderson, W. Bates, F. Salmassi, P. Nachimuthu, A. Pearson, D. Richardson, D. Hambach, N. Hoffmann, W. Hasse, and K. Hoffmann, "Electromigration in integrated circuit interconnects studied by X-ray microscopy," *Nuclear Instruments & Methods in Physics Research Section B*, **199**, pp. 469–474, 2003.
15. ALS Highlights 2003, available at <http://www-als.lbl.gov/>
16. D. Weiß, G. Schneider, B. Niemann, P. Guttmann, D. Rudolph and G. Schmahl, *Ultramicroscopy* **84**, (2000) 185–197.
17. G. Schneider, *Ultramicroscopy* **75**, (1998) 85–104.
18. S. Vogt et al., *Journal of Structural Biology* **132**, (2000) 123–132.
19. W. Meyer-Ilse et al., *Journal of Microscopy* **201**(3), (2001) 395–403.
20. G. Schneider et al., *Surf. Rev. Lett.*, Vol. 9, No. 1 (2002) 177–183.
21. Larabell, C. A. and Le Gros, M. A. (2003). X-ray Tomography Generates 3-D Reconstructions of the Yeast, *Saccharomyces cerevisiae*, at 60 nm Resolution. *Molecular Biology of the Cell*. Submitted.

## **PART 2**

# **Further Development of Soft X-ray Scanning Microscopy with an Elliptical Undulator at the Advanced Light Source**

---

TONY WARWICK,<sup>1</sup> HARALD ADE,<sup>2</sup> SIRINE FAKRA,<sup>1</sup> MARY GILLES,<sup>1</sup> ADAM HITCHCOCK,<sup>3</sup> DAVID KILCOYNE,<sup>2</sup> DAVID SHUH,<sup>1</sup> AND TOLEK TYLISZCZAK<sup>1</sup>

<sup>1</sup>Lawrence Berkeley National Laboratory, Berkeley, CA 94720, USA

<sup>2</sup>Department of Physics, North Carolina State University, Raleigh, NC 27695, USA

<sup>3</sup>Brockhouse Inst. Mtrls. Res., McMaster University, Hamilton, ON L8S 4M1, Canada

Soft X-ray scanning microscopy (1) is under continuing development at the Advanced Light Source. Significant progress has been made implementing new scan control systems in both operational microscopes (2) and they now operate at beam lines 5.3.2 and 11.0.2 with interferometer servo scanning and stabilization. The interferometer servo loop registers the images on a universal x/y coordinate system and locks the X-ray spot on selected features for spectro-microscopic studies. At the present time zone plates are in use with 35 nm outer zone width and the imaging spatial resolution is at the diffraction limit of these lenses.

Current research programs are underway in areas of polymer chemistry, environmental chemistry and materials science.

A dedicated polymer STXM is in operation on a bend magnet beam line (4) and is the subject of a separate article (3) in this issue. Here we focus on the capabilities of STXM at a new beam line that employs an elliptical undulator (5) to give control of the polarization of the X-ray beam. This facility is in the process of commissioning and some results are available, other capabilities will be developed during the first half of 2003.



Figure 1: Beamline 11.0.2 at the ALS. The monochromator is foil-wrapped at the lower left. The line branches by means of a movable mirror and each branch has its own set of exit slits. The STXM branch is seen here, with the STXM enclosure to the upper right.

### Beamline capabilities

Figure 1 shows the beam-line. A newly engineered entrance-slit-less SX700 monochromator operates in collimated light (6,7) and serves the microscope part-time. The new monochromator has been engineered to emphasize cooling and cleanliness. The flux reduction at the carbon edge is about 30%. Much of the research on this line will involve organic molecules, studied at the C1s absorption edge.

Figure 2 is the standard resolution test, measuring  $N_2$  vibrational levels using a partial pressure of air in the microscope enclosure. Figure 3 shows the computed count rates in the microscope, operating the monochromator with three different values of the focusing parameter  $C_{ff} = \cos \beta / \cos \alpha$ .

$C_{ff} = 1.3$  has higher diffraction efficiency and smaller dispersion for maximum count rate.  $C_{ff} = 3.9$  provides high spectral resolution with a large safety margin against vertical electron beam motion.

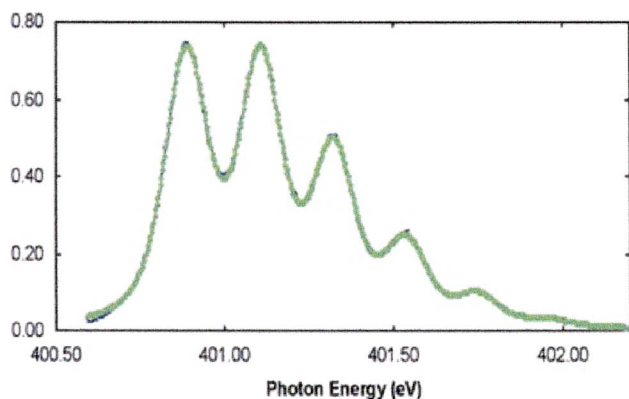


Figure 2:  $N_2$  absorption spectrum, with a spectral resolving power of around 8000.

Microscope Flux at  $R=7500$ , 150 lines/mm (photons/second)

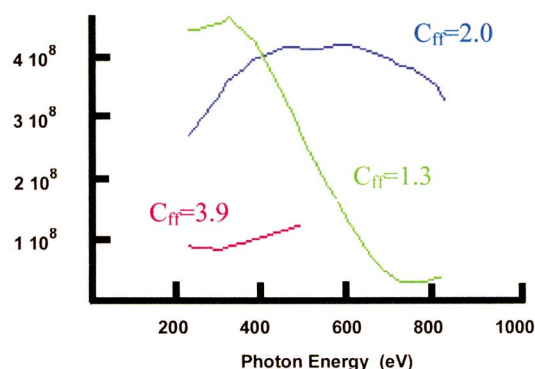


Figure 3: Computed microscope illumination.

### Microscope performance

At this time the microscope operates at room temperature in an enclosure that can be evacuated, or operated with a helium atmosphere. X-rays are focused by a zone plate with 155  $\mu m$  diameter with 35 nm outer zones (8). The instrument operates with photon energy from 150 eV to 2000 eV. The detector is a phosphor-photomultiplier combination detecting transmitted photons with an efficiency between 5% and 50%. Figure 4 shows the spatial resolution achieved with the microscope at the new beam line. In addition to the advantages of the interferometer servo system the new beam line offers higher count rates and lower intensity fluctuations, so that better quality images can be obtained at faster acquisition rates.

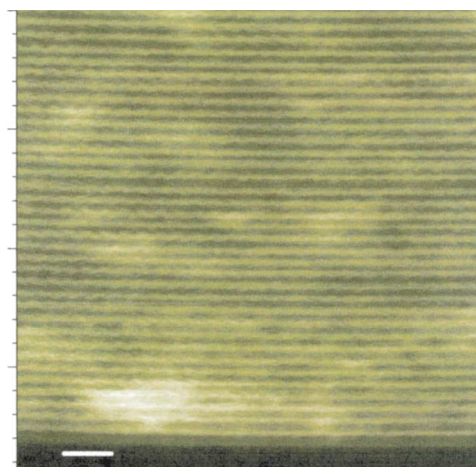


Figure 4: X-ray image of the outer zones of a zone plate with 25 nm outer zone width. The bar in the figure is 200 nm long. This image was acquired counting for 1 ms per pixel at 390 eV.

#17

Physical Sciences Laboratory

1/2 pg. Horiz.

Repeat SRN 16.2

---

## Linear Dichroism Microscopy of Semi-Crystalline Polymers

D. KILCOYNE (NCSU) T. TYLISZCZAK (LBNL), H. ADE (NCSU),  
Y. WANG, AND M. RAFAILOVICH (SUNY@STONY BROOK)

Semi-crystalline polymers, including most thermoplastics, are a class of materials that are ubiquitous, yet various aspects of their properties, including their crystallization behavior, are still poorly understood. In thin films, due to the effects of the interfaces, the issues are compounded to the point that it is unclear to what degree the films crystallize and what orientation crystallites have. X-ray linear dichroism microscopy (9) is a potentially very powerful tool to address many open questions, but dichroism microscopy applications have been limited due to difficulties with image registration if samples are rotated. The elliptical undulator that illuminates the 11.0.2 STXM will eliminate these constraints. Complete control over the orientation of linear polarization will allow the characterization of the orientation of molecular bonds. First results with controlled linear polarization indicate that the 11.0.2 STXM will become a powerful tool for semi-crystalline polymeric materials.

Figure 5 shows spherulitic-like domains in a thin poly-ethylene film. The relative contrast of these domains changes as the linear polarization of the EPU is switched from horizontal to vertical. Since the sample and microscope are completely unchanged, complete registration between images that contain orientation information can be achieved. A detailed analysis of images such as these and complete NEXAFS spectra will allow the determination of important parameters of the sample.

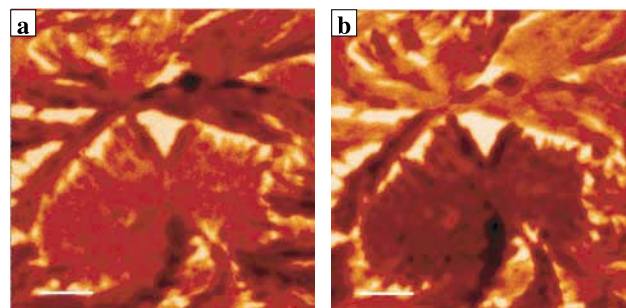


Figure 5: Thin poly-ethylene film with spherulitic crystallites imaged at 288 eV with horizontal (left) and vertical (right) linearly polarized X-rays. (Scale bar = 5  $\mu$ m).

## Molecular Environmental Science: Metal Bio-accumulation in Bio-films

T. TYLISZCZAK (LBNL), A. P. HITCHCOCK AND  
T. ARAKI (McMASTER UNIVERSITY), J. R. LAWRENCE  
(NWRI, SASKATOON), AND G. G. LEPPARD (NWRI, BURLINGTON)

Bio-films consist of colonies of microbial cells that use exo-cellular polysaccharides, proteins and nucleic acids to form an optimized environment. We are studying biofilms in order to understand how they deal with toxic elements such as metals or organic species they encounter in their environment. In addition to fundamental understanding, we hope this research will lead to practical strategies for biofilm-based remediation of polluted sites. The new undulator-based 11.0.2 STXM is an ideal tool for such studies since it can be applied to fully wet samples, it covers a wide spectral range, and it has high sensitivity on account of excellent energy resolution and high brightness.

Figure 6 shows our first results from measurements with STXM 11.0.2. Here various metals in a river water biofilm are mapped relative to the biochemical structure of the cells and surrounding exo-cellular matrix. The Ni has been localized by bacteria from a 24-hour exposure to a flowing stream of 1 ppm nickel chloride.

## Materials Science: Studies for the Formation of Silicon Nanocrystals

D. GERION AND N. ZAITSEVA (LLNL),  
S. FAKRA (LBNL)

The wet chemistry synthesis of colloidal nanocrystals with tunable size and shapes offers exciting possibilities in material science. Nanocrystals are toolkits for the fabrication of nanodevices and nanosensors. They also allow the study of properties related to quantum size effects. Well-established examples are the gram quantity synthesis of II–VI semiconductor quantum dots (CdSe, CdTe), of magnetic Co dots and disks, and of ferroelectric oxides. Silicon nanocrystals would be technologically important, but have not yet been produced on a large scale.

Following an idea of Korgel et al. (10), we investigated the thermal decomposition of organo-silane precursors at high temperature and high pressure. As an example, we thermally decomposed diphenylsilane (Phenyl<sub>2</sub>-Si-H<sub>2</sub>) in hexane at 400°C and >100 bars. Due to these extreme conditions, extensive degradation of both precursor and solvent occurs. This generates new chemical species that we wish to analyze at the molecular level.

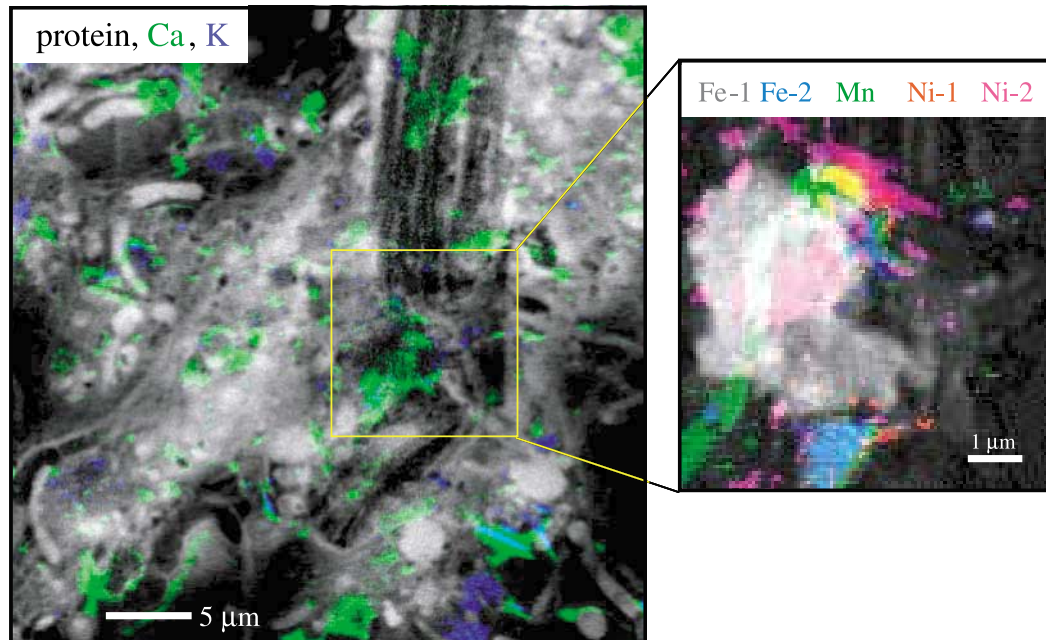


Figure 6: STXM-derived chemical maps of a river water biofilm. (left) Grayscale image of biological components (cells, exocellular matrix) determined from (I288-I279), the difference in images recorded at 288.2 eV (amide  $p^*$  band) and 279 eV (pre-C1s), overlaid with maps of  $\text{CaCO}_3$  (I352-I350) and  $\text{K}^+$  (I294-I293). (right) Detailed study of metal distributions in the indicated region. Image sequences recorded in the metal 2p spectral regions show multiple metal environments for Fe, and Ni. The color coded composite maps these components as well as Mn.

#1

ACCEL Instruments

1/2 pg. Horiz.

Repeat SRN 16.2

Figure 7 shows the Si K-edge absorption spectra of a pure Si compound and of the result of the thermal decomposition of diphenylsilane at 400°C. The shift in the absorption edge indicates a different bonding environment for Si atoms. Tests on some molecular precursors (organo-silane compounds) at the C, O and Si K absorption edges using BL11.0.2 STXM are underway to see whether we can break the Si-C bonds and cause the condensation of Si ions into pure Si nanoparticles.

#### Future directions

The microscope installed at beam line 11.0.2 covers the C1s edge, where much work is done. It can also work at higher energy. The transition metal L edges are important environmental contaminants whose heterogeneous chemistry is of great interest in a wide variety of media, particularly with water present.

The availability of elliptically polarized radiation will spawn a new program in magnetic microscopy, not yet underway. The M edges of rare earth magnetic materials will be accessible. The higher flux and high photon energy will make fluorescence studies a possibility for future development. ■

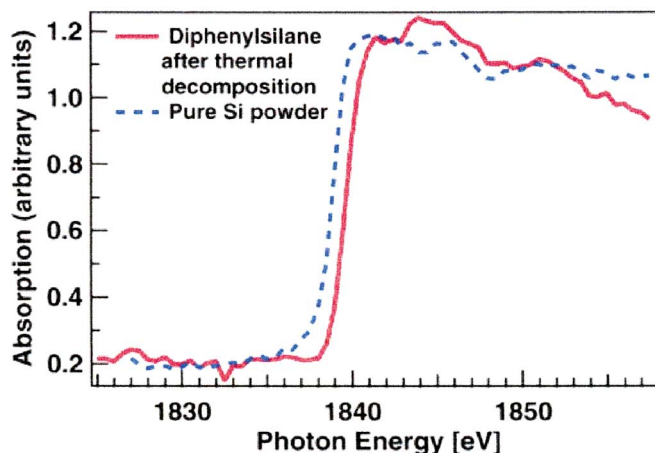


Figure 7: NEXAFS spectra extracted from image stacks on two test compounds, acquired with the 11.0.2 STXM. These are silicon K-edge absorption spectra of a pure Si compound and of the result of the thermal decomposition of diphenylsilane at 400°C. Diphenylsilane is composed of two phenyl rings and two hydrogen atoms attached to a central Si atom. These chemical shifts will allow high contrast images of heterogeneous silicon chemistry.

# Acknowledgements

Zone plates used at 11.0.2 are provided by the Center for X-ray Optics, LBNL, in a collaboration with B. A. Harteneck, D. L. Olynick, J. A. Liddle, and E. Anderson. The research and development efforts of this group are essential to the scanning microscopy programs at the ALS. We also thank all the members of the ALS vacuum group, under the supervision of D. Colomb, for their effort and expertise.

# References

1. J. Kirz and H. Rarback, "Soft X-ray Microscopes," *Rev. Sci. Instrum.* **56**, 1 (1985).
2. A. L. D. Kilcoyne, T. Tyliczszak, W. F. Steele, S. Fakra, P. Hitchcock, K. Franck, E. Anderson, B. Harteneck, E. G. Rightor, G. E. Mitchell, A. P. Hitchcock, L. Yang, T. Warwick and H. Ade, "Interferometer Controlled Scanning Transmission X-ray Microscopes at the Advanced Light Source," *J. Synchrotron Rad.* **10**, 125 (2003).
3. H. Ade, A. P. Hitchcock, G. E. Mitchell, A. L. D. Kilcoyne, T. Tyliczszak, R. Fink and T. Warwick, "Applications with the Dedicated Scanning Transmission X-ray Microscope at ALS Beamline 5.3.2," *Synchrotron Radiation News*, this issue.
4. T. Warwick, H. Ade, D. Kilcoyne, M. Kritscher, T. Tyliczszak, S. Fakra, A. Hitchcock, P. Hitchcock and H. Padmore, "A New Bend Magnet Beam Line for Scanning Transmission Microscopy at the Advanced Light Source," *J. Synchrotron Rad.* **9**, 254 (2002).
5. S. Marks, C. Cortopassi, J. DeVries, E. Hoyer, R. Leinbach, Y. Minamihara, H. Padmore, P. Pipersky, D. Plate, R. Schleuter and A. Young, "The Advanced Light Source Elliptically Polarizing Undulator," *Proc. Particle Accel. Conf.*, Vancouver, Canada, (1997).
6. R. Follath and F. Senf, "New Plane Grating Monochromators for Third Generation Synchrotron Radiation Light Sources," *Nucl. Inst. Meth. Phys. Res.* **A390** 388 (1997).
7. T. Warwick, D. Shuh, "MES project (phase 1) conceptual design," ALS internal note LSBL 536, October 1999.
8. E. H. Anderson, D. L. Olynick, B. Harteneck, E. Veklerov, G. Denbeaux, W. L. Chao, A. Lucero, L. Johnson, and D. Attwood, "Nanofabrication and Diffractive Optics for High-Resolution X-ray Applications," *J. Vac. Sci. & Technol. B* **18**, 2970 (2000).
9. H. Ade and B. Hsiao, "X-ray Linear Dichroism Microscopy," *Science* **262**, 1427 (1993).
10. J. D. Holmes, K. J. Ziegler, R. C. Doty, L. E. Pell, K. P. Johnston and B. A. Korgel, "Highly Luminescent Silicon Nanocrystals with Discrete Optical Transitions," *J. Am. Chem. Soc.*, **123** 3743 (2001).

#19

Taylor & Francis  
House Ad

1/2 pg. Horiz.

NEW